sonic

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sonic  
Dual action  
sonic technology

## INTRODUCTION

**DENTAID** is an international company and pioneer in research, development and marketing of dental care products. Founded with the aim of improving dental health worldwide, we offer the best solutions endorsed by both dental and pharmaceutical professionals.

**VITIS** is the DENTAID brand for **daily oral hygiene** products to help in **prevention, care and effective protection of your mouth**. In order to ensure good oral health, we have gradually developed **the most comprehensive and specialised range** of toothbrushes, toothpastes, mouthwashes and dental flosses and tapes, adapted to your needs in every case and at any time.

The **VITIS brand specialises in the development of high quality manual and electric toothbrushes**. It features a wide range of brushes designed by experts in different areas to cover the different needs of your mouth, providing effective removal of dental biofilm, while caring for gums, teeth and other oral tissues.

Not only are **VITIS sonic S10 & S20** electric toothbrushes backed by the experience of a specialised toothbrush brand, but they also offer cutting-edge sonic **technology for superior brushing**.

## THE IMPORTANCE OF GOOD TOOTHBRUSHING

Over 700 bacterial species live in our mouth, and can grow and multiply on the different oral tissues. When these bacteria adhere to the tooth's surface, they form a film called **dental biofilm**. This **biofilm is responsible for the appearance of different diseases and oral disorders such as tooth decay, periodontal diseases including gingivitis and periodontitis, and halitosis**.

**Regular control of oral biofilm is therefore essential to maintain proper oral health**. A toothbrush is the basis for good daily oral hygiene. A quality brush used with a proper brushing technique applied for the right amount of time, will remove the biofilm and food debris deposited on the teeth accomplishing proper cleanliness of teeth and gums.

**Electric toothbrushes with sonic technology have been shown to enhance oral biofilm removal by improving brushing technique, brushing time compliance and increased user motivation for thorough oral hygiene**.

## 1. VITIS SONIC S10 / S20 TECHNOLOGY FOR EFFECTIVE BRUSHING

The VITIS sonic electric toothbrush, with sonic technology, emits high frequency acoustic vibrations which cause lateral oscillation of the filaments at a speed of up to 40,000 oscillations per minute.

The basic advantage of sonic brushing is that it adds the **benefits of the dual action of sonic technology** to the manual action of mechanical cleaning:

- **Mechanical action:**

the high frequency vibration of the filaments on the surface of the teeth favours the displacement of the dental biofilm leading to its removal.

- **Hydrodynamic action:**


- The high speed vibration of the filaments emits acoustic waves; these waves and their energy are transmitted through the saliva and throughout the oral cavity, attaining greater dissolution of the oral biofilm.
- The vibration generates small bubbles of oxygen discharged into the saliva and which displace the bacterial plaque in difficult to reach areas.

VITIS sonic **heads have a unique design** that boosts the benefits of brushing, providing greater precision in removing the biofilm, attaining difficult to reach areas such as the interdental zones and the gingival margin, while nurturing both teeth and gums.

All thanks to the following **characteristics of the VITIS heads:**

- **Small size or access:** favours cleaning hard-to-reach areas, interdental spaces and the gingival margin.
- **High quality Tynex® filaments** with rounded tips, for efficient brushing, while protecting teeth and gums.
- **Exclusive design and distribution of the filaments** in the head means we offer heads adapted to the needs of each user.
- **Monotip zone:** clearing even the most resistant dental biofilm.
- **A protective cap** to shield the filaments from contact with the exterior but allowing the passing of air for quick drying.

**Thanks to their sonic technology and exclusive design, VITIS sonic toothbrushes are an ideal solution** as they eliminate plaque efficiently, even in hard-to-reach places; they manage to improve brushing technique and support compliance with brushing time, thereby protecting both teeth and gums.

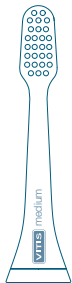**VITIS sonic** is recommended for **everyone**, and especially for those who:

- Demand high standards of dental hygiene
- Need to improve their brushing technique
- Have plaque build-up
- Have problems with their gums
- Are undergoing orthodontic treatment
- Have dental implants
- Have limited manual dexterity (arthrosis, arthritis, etc)

## 2. DESCRIPTION OF THE VITIS SONIC ELECTRIC TOOTHBRUSH

The VITIS sonic toothbrush range includes models S10 & S20 offering sonic technology with different features.

### **VITIS sonic S10**


Maximum brushing efficacy while protecting teeth and gums

- A** **Brush head**, high quality Tynex® filaments with rounded tips for efficient brushing and nurturing of teeth and gums
- B** **Protective cap**; keeps the filaments together and in peak hygienic condition
- C** **Coloured ring**; allows differentiation of the head if there are several users
- D** **Ergonomic handle** with non-slip grooves
- E** **On / Off Button**
- F** **Timer** with quadrant change indicator (30 seconds per quadrant)
- G** **Charge indicator** light indicates whether the toothbrush is properly charged
- H** **Charging cradle** compatible with both VITIS sonic models


## VITIS sonic S10

Maximum brushing efficacy  
while protecting teeth and gums
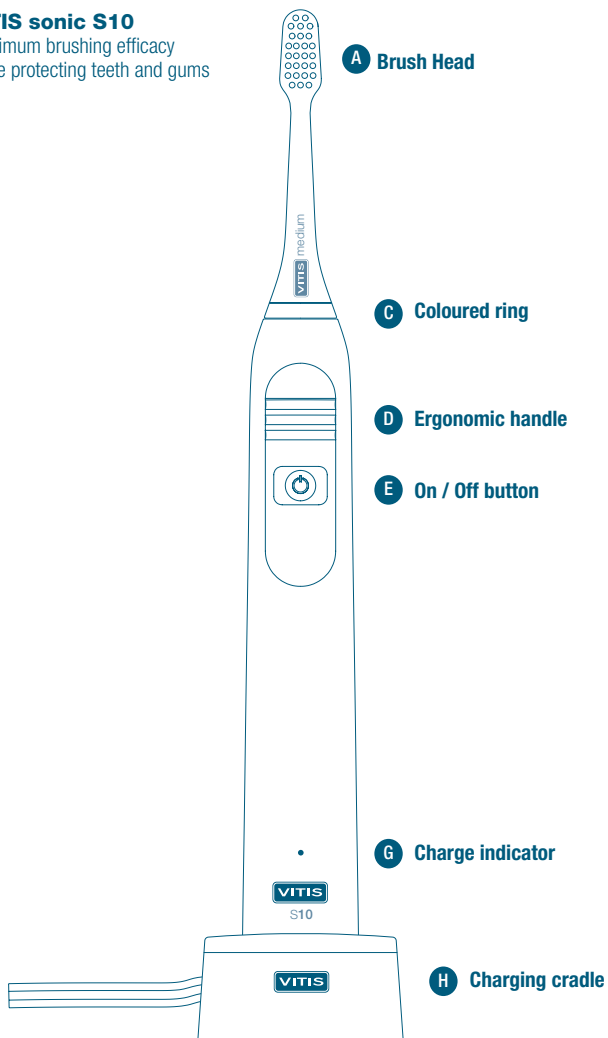
## VITIS sonic S20

Maximum precision and efficiency for superior brushing while protecting teeth and gums

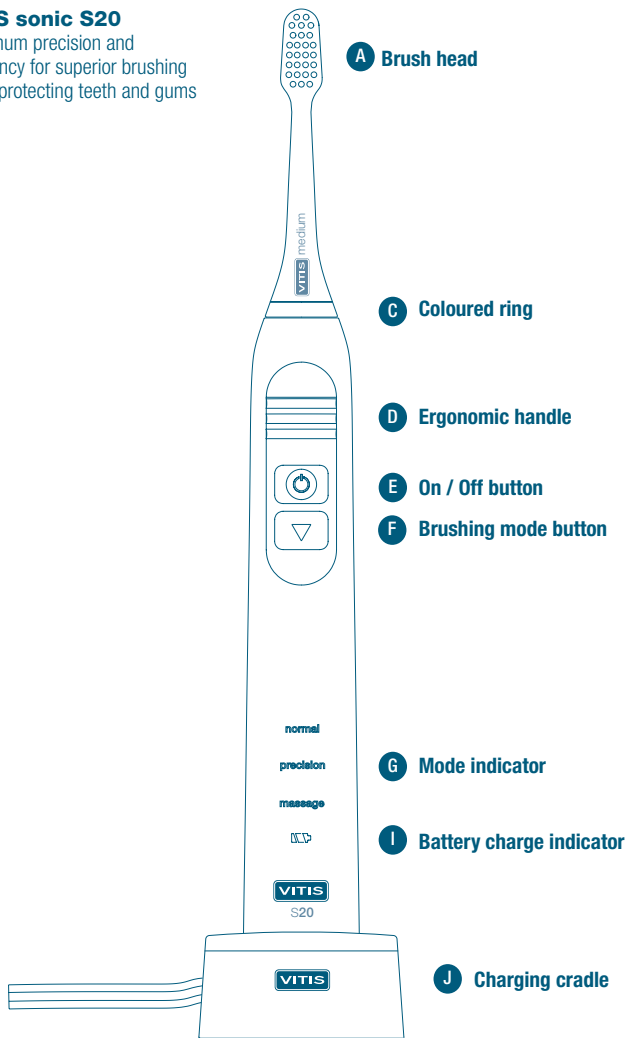- A** **Brush head**, high quality Tynex® filaments with rounded tips for efficient brushing and nurturing of teeth and gums
- B** **Protective cap**; keeps the filaments together and in peak hygienic condition.
- C** **Coloured ring**; allows differentiation of the head if there are several users
- D** **Ergonomic handle** with non-slip grooves
- E** **On / Off Button**
- F** **Brushing mode button: Normal · Precision · Massage**
- G** **Mode indicator**, displaying which brushing mode is in use
- H** **Timer** with quadrant change indicator (30 seconds per quadrant)
- I** **Battery charge indicator**
- J** **Charging cradle** compatible with both VITIS sonic models

### 3. GETTING STARTED



When using the VITIS sonic toothbrush for the first time, perform the following steps:

1. Place the brush head on the handle, with the filaments facing forward, and press down until it locks.
2. Plug the charger into the wall socket.
3. Place the toothbrush on the charging cradle.

**S10:** The light shows (charge indicator) that the toothbrush is correctly placed and is charging.

**S20:** The light (battery level indicator) shows the charge level of the electric toothbrush battery.

4. The battery will be fully charged after being plugged in for 24 hours in the cradle. It is recommendable the battery be left to charge for 24 hours if:
  - It is the first time it is being used
  - It has been left unused for over three months
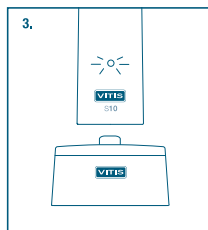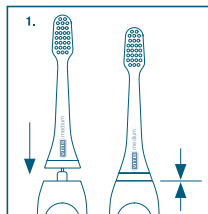
### 4. RECOMMENDED USE

To achieve maximum efficacy of your VITIS Sonic brush, follow the steps that appear below. This will also allow you not only to remove plaque from teeth (dental biofilm), but also to remove plaque from the gingival margin between teeth and gums, which causes the onset of gum diseases such as gingivitis and periodontitis.

1. Apply toothpaste to the brush filaments.
2. Insert the toothbrush in your mouth placing the filaments along the teeth at a 45-degree angle to the gumline.
3. Press the on/off button.
4. When the VITIS sonic brush vibrates gently against the gumline, sweep the filaments away from the gums.

- Brush the outer and inner surfaces of all teeth in this manner.
- Clean chewing, or occlusal, surfaces using horizontal movements.
- Brushing should last at least 2 minutes: 30 seconds per quadrant. A brief pause in the brush movement will indicate when to move on to the next quadrant. Once brushing has continued for 2 minutes, the brush will completely stop.

- If you wish to change the brushing mode (only possible for the VITIS Sonic S20), press the 'change mode' button. The brush handle displays which brushing mode you are using.

5. When the VITIS Sonic brush stops vibrating, the recommended brushing time has finished. If you wish to continue using the brush, press the on/off button again.

## 5. AVAILABLE MODES (ONLY VITIS SONIC S20 MODEL)

The **VITIS sonic S20** brush offers 3 brushing modes to choose from, depending on the user's preferences or needs:

- **Normal:** 31.000 strokes per minute, allowing for maximum efficacy in the removal of oral biofilm for proper daily oral hygiene and care of teeth and gums.
- **Precision:** 40.000 strokes per minute. More strokes per minute make it so that the filaments have less of a stretch, resulting in greater precision and maximum efficacy in the removal of oral biofilm, achieving deep cleaning in a more localised area.
- **Massage:** Combined speeds (31,000 and 40,000 strokes/min) to clean and massage the line that separates gums and teeth and to gently remove plaque from all oral tissues.

## 6. TIMER

Both VITIS Sonic brush models have a timer that pauses brushing briefly every 30 seconds to indicate when to move on to the next quadrant.

After operating for the recommended brushing time period of two minutes, the brush stops permanently, indicating brushing completion. This way you can be sure you have fulfilled the recommended brushing time for optimum oral cleaning of all sides of all teeth.

## 7. AVAILABLE BRUSH HEADS

The VITIS sonic S10 and S20 electric brushes are compatible with the following interchangeable brush heads:

- **VITIS sonic medio/medium brush head:** with medium strength filaments for maximum removal of bacterial plaque while gentle on teeth and gums.
- **VITIS sonic encías/gingival brush head:** with extra-soft filaments. Intended for people with delicate gums, with signs of inflammation or people who prefer softer brushing. Allows for deep cleaning without damaging gums and with an extremely smooth feel.

## 8. MAINTENANCE

For proper upkeep of the brush and its head, it is advisable to:

1. Turn off the brush and separate the head from the handle after use.
2. Rinse the brush head with water and dry it.
3. Clean the handle and dry it too. Once it is dry, store it on the charging cradle. Clean the charging cradle regularly with a damp cloth, first making sure that the charging cradle is unplugged.

## 9. SAFETY AND WARNINGS

The VITIS sonic S10 and S20 electric toothbrushes are safe and effective, provided they are used according to the instructions. Like any electrical appliance, you should follow the safety precautions, especially when children are present.

Before using the product read the instructions and keep this manual for future reference.

To reduce the risk of burns, electrocution, fire or injury to any person:

### **WATER – danger**

- Do not immerse the charger in water or other liquids.
- Do not place the charger on or near a bathtub or sink with water, where it can fall.
- Do not attempt to retrieve the charger if it has dropped into water. Unplug it immediately.
- Do not use the unit while bathing.
- Unplug the charger from the wall socket to clean it and make sure that it is completely dry before connecting it again. Cleaning and maintenance should not be performed by children.

### **VOLTAGE – danger**

- Before turning the unit on, check to make sure the voltage specification (shown on the charging cradle) matches the local mains voltage.
- This product is designed to operate in a range of 100-240 volts


- Never force the plug into an outlet.
- If the plug does not fit properly into the socket, do not use it.

## **SAFETY WARNINGS**

- This product is not intended for use by persons (including children) with limited physical, sensory or mental capabilities or with a lack of experience or knowledge, unless they have been properly trained to use the electric toothbrush safely.
- Keep out of the reach of children if unsupervised by an adult. This product is not designed to be used by children under the age of 3 years. Children should not play with the product.
- Children between the ages of 3 and 14 years and people with physical, sensory or mental disabilities or who lack experience or knowledge, may use the brush while supervised or if the instructions for use of the brush have been carefully explained and the risks involved have been understood.
- Use the product only as described in this manual. Do not use accessories or brush heads other than those recommended by the manufacturer.
- Do not use the product if the power cord, the base or the charger are damaged. Do not use the product if it is not working properly. Stop using the product immediately if you notice that it is damaged. Do not tamper with or attempt to repair the product; this may cause fire, electric shock or injury.
- If the supply cord is damaged, it must be replaced by the manufacturer, its service agent, or similarly qualified persons in order to avoid a hazard.
- Keep the cord away from sources of extreme heat.
- Do not put the brush in the dishwasher.
- When the charger is connected to the mains, keep it away from hot surfaces and out of direct sunlight.
- The battery must be removed from the appliance before it is scrapped. The appliance must be disconnected from the mains supply when removing the battery. The battery is to be disposed of safely. The consumer must use the appropriate tools to remove the screws from the product, then disengage the clip fittings on the outer casing to access the battery for removal, prior to the product being scrapped.

## 10. END-OF-LIFE DISPOSAL (ENVIRONMENT)

The brush handle contains a recyclable Li-ion battery. To preserve the environment, at the end-of-life of the product, remove the battery before disposing of the brush, or take it to a collection point provided for that purpose.
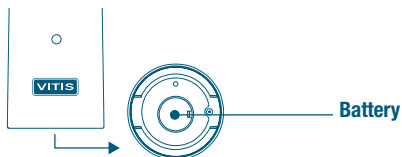
## 11. WARRANTY

DENTAID offers a two-year warranty from the date of product purchase. Should you detect that the brush has a material- or manufacture-related defect within the warranty period, DENTAID promises to replace it with a new brush. The warranty does not cover defects deriving from misuse, tampering or damage after purchase.

Proof of purchase of the product showing the date the product was obtained must be submitted in order for this warranty to be considered valid.

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**LABORATORIOS DENTAID S.A.**